CODE No. 4700

MBL

Released: July 31, 1997

Revised: January 10, 2012

For Research Use Only. Not for use in diagnostic procedures.

For the sensitive detection and quantification of early-stage apoptosis

MEBCYTO[®]

Apoptosis Kit (Annexin V-FITC Kit)

CODE No. 4700 100 tests

MBL MEDICAL & BIOLOGICAL LABORATORIES CO., LTD.

URL <https://ruo.mbl.co.jp>

e-mail support@mbi.co.jp, TEL 052-238-1904

CONTENTS

1. Intended Use -----	1
2. Summary and Explanation-----	1
3. Principle -----	1
4. Materials provided -----	2
5. Storage and Stability -----	2
6. Materials and equipment required -----	2
7. Precautions -----	2
8. Procedure -----	2
9. Related products-----	4
10. References -----	4

Before use, thoroughly read these Instructions.

Intended Use

The MEBCYTO[®] Apoptosis Kit (Annexin V-FITC Kit) is intended for the detection of early stage of apoptotic cells in research samples. **For Research Use Only. Not for use in diagnostic procedures.**

Summary and Explanation

Apoptosis, a term that describes regulated cell death, is fundamental feature of many biological processes including normal development, homeostasis, and disease. Early during the process of apoptosis, cells lose their phospholipids membrane asymmetry and expose phosphatidylserine (PS) at the cell surface. This process can be monitored by using Annexin V which is a Ca²⁺ dependent, phospholipids-binding protein with high affinity for PS, and is useful for identifying apoptotic cells with exposed PS.

Translocation of PS to the external cell surface is not unique to apoptosis, but occurs also during cell necrosis. The difference between these two forms of cell death is that during the initial stage of apoptosis the cell membrane remains intact, while at the very moment that necrosis occurs the cell membrane loses its integrity and becomes leaky. Therefore, necrotic cells easily stained with Propidium Iodide (PI) as well as Annexin V, where as Apoptotic cells stained only with Annexin V.

The MEBCYTO[®] Apoptosis Kit (Annexin V-FITC Kit) can identify cells in an earlier stage of apoptosis than assays based on DNA fragmentation.

Principle

When staining the cells with Annexin V-FITC and Propidium Iodide (PI), it can be used in a bivariate analysis to distinguish between cells undergoing apoptosis (PI negative) and those that are necrotic or dead (PI positive).

Materials provided

Each kit contains:

Name	Materials	Quantity (for 1 kit)
Annexin V-FITC	FITC conjugated Annexin V (ready-for-use)	1.0 mL × 1 vial
Propidium Iodide (PI)	Propidium Iodide (PI) (100 µg/mL) (ready-for-use)	0.5 mL × 1 vial
Binding Buffer	Binding buffer (ready-for-use)	50 mL × 1 bottle

These reagents are sufficient for 100 tests.

Storage and Stability

All kit components must be stored at 2-8°C. All reagents are stable for 24 months after manufacturing when stored at the indicated conditions.

Materials and equipment required

- Phosphate buffered saline (PBS)
- Adequate test tube
- Adjustable micropipette
- Fluorescence activated cell sorter

Precautions

1. **Annexin V-FITC** contains sodium azide (0.09%) as preservative. Azide may react with copper or lead in plumbing system to form explosive metal azide. Therefore, always flush with plenty of water into a drain when disposing materials containing azide.
2. **Propidium Iodide (PI)** is highly toxic. Handle with care.
3. This kit is intended for research use only. Not for use in diagnostic procedures.

Procedure**◆ Preparation of Reagents**

All reagents are ready for use.

◆ Assay procedure

Since apoptosis is a very rapid and dynamic process, we recommend that you perform analysis immediately after staining.

STEP A. (Sample Staining)

- 1) Induce apoptosis by desired method.
- 2) Wash the cells (2.0×10^5 cells) to be used for staining once with PBS.

Note

For adherent cells, trypsinize and gently wash cells once with serum-containing medium followed by wash once with PBS.

- 3) Resuspend the cells in 85 μ L of **Binding Buffer**.
- 4) Add 10 μ L of **Annexin V-FITC** and 5 μ L of **Propidium Iodide** to 85 μ L cell suspension prepared as given by **STEP A.-3**.
- 5) Mix well and incubate at room temperature (20-25°C) for 15 minutes in the dark.

STEP B. (Analyzing the stained cells by Flow Cytometry)

- 1) After incubation (**STEP A.-5**), add 400 μ L of **Binding Buffer**.
- 2) Measure the cell sample by Flow Cytometry using a single laser emitting excitation light at 488 nm.

Note

The signal generated by Annexin V-FITC can be detected in the FITC signal detector (usually FL1). The signal generated by Propidium Iodide can be monitored by the detector for phycoerythrin emissions (usually FL2 or FL3). In case the multicolor assay, careful calibration of the flow cytometer is needed to compensate for the overlapping emission.

◆ Example of assay result

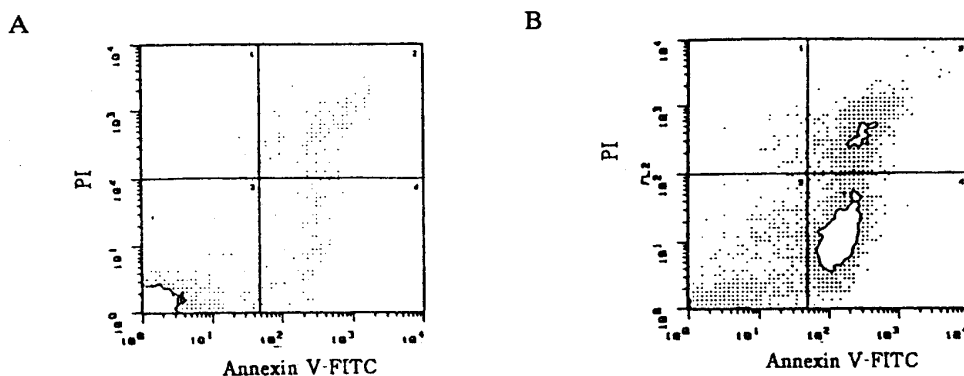
Fig. 1 Example of assay result

A: Viable Jurkat cells (negative control).

B: Flow Cytometry analysis of the Jurkat cells induced apoptosis by anti-Fas antibody (MBL; code no. SY-001).

CODE No. 4700

Related products

4690	APOPCYTO Annexin V-Azami-Green Apoptosis Detection Kit
8441	MEBSTAIN Apoptosis TUNEL Kit II
8445	MEBSTAIN Apoptosis TUNEL Kit Direct
SY-001	anti-human Fas (CD95) (CH-11)

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